

Please download Fiji at:

fiji.sc

and demo data, scripts, slides at:

<https://github.com/ABIF-McGill/ABIF-Fiji-Workshop-2025>

Fiji workshop 2025

- Day 1:
 - Image data visualization in Fiji
- Day 2:
 - Quantitative image analysis in Fiji
- Day 3:
 - One-on-one analysis support with your own image data

Joel Ryan, Philip Kesner, Barbara Da Rocha, Mar Garcia Ferrés


ABIF

March 4 to 6 2025

Day 1: Data visualization in Fiji

- Today's goals
 - Get acquainted with Fiji!
 - LUTs, contrast
 - Composite images
 - Intro to scripting!
 - Look at data from 3D images
 - Look at data from timelapse images
 - Images for publications


Fiji is awesome! But it has some drawbacks...

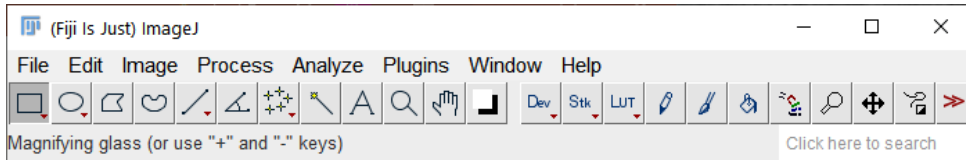

- Why should you use Fiji?
 - Free open-source software
 - Huge development community
 - Libraries to open nearly any microscope image (!)
 - Automation and scripting with multiple programming languages (!)
- When should you use something else than Fiji?
 - More powerful computation
 - Better 3D visualization
 - Many deep learning tools
 - Instrument-specific processing

Digital images

- A digital image is simply a table of numbers
- Each number represents the intensity value of a pixel
- The intensity value is proportional to the amount of light collected at that location
- Software such as Fiji simply displays each number as a certain shade or colour

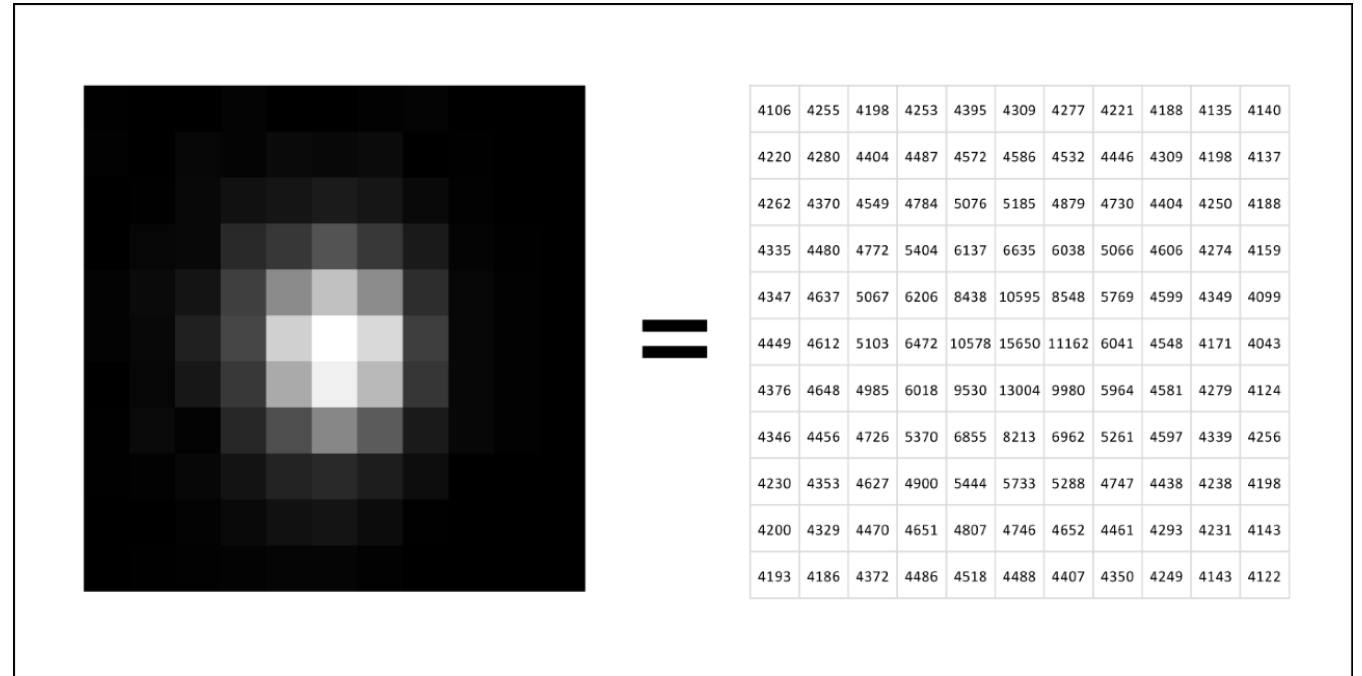

Example images


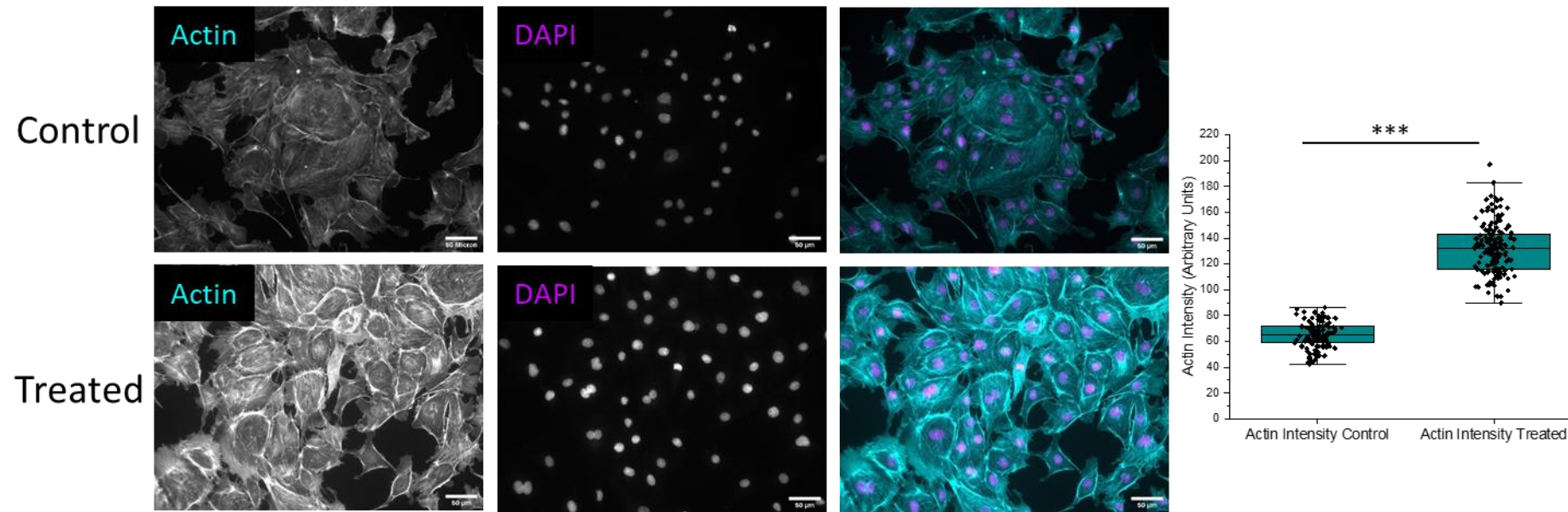





Figure 1: Images of CHO-K1 cells stained with Phalloidin Alexa Fluorophore 488 and stained with the nuclear probe DAPI. Images were flat field corrected using images of a fluorescent plastic slide imaged on the same microscope. Images were collected on an upright Zeiss Axioskop with a EC PlanNeoFluar 20x/0.5 NA objective lens with an AxioCam ICm1 camera. DAPI was imaged using a DAPI cube and an exposure time of 50 ms. Actin was imaged with a FITC cube and a 150 ms exposure. Brightness, gamma and contrast were adjusted to visualize the actin features well. The scale bar is 50 µm. Box plot of actin intensity for control (n=110) and treated (n=149) cells. T-test was done with a two sample t-test, unequal variance, two tailed. *** corresponds to $P < 0.001$.

Example images




Imaris demo image

Drosophila egg Chamber
3 channels x 98 z-slices
16 frames-per-second playback

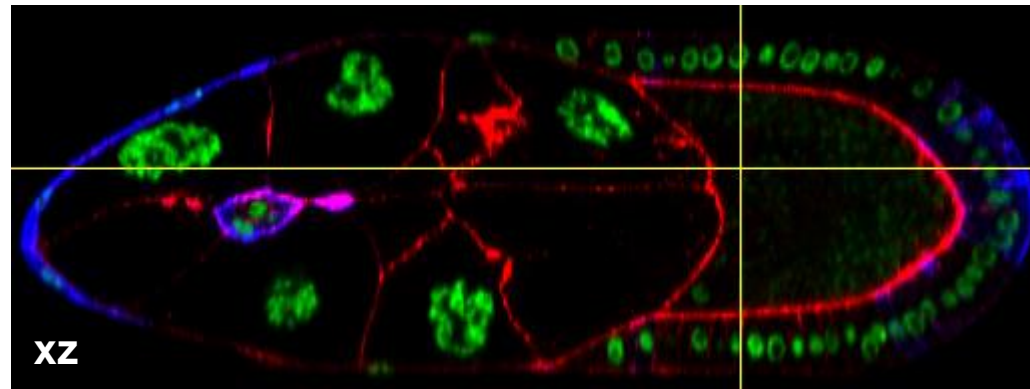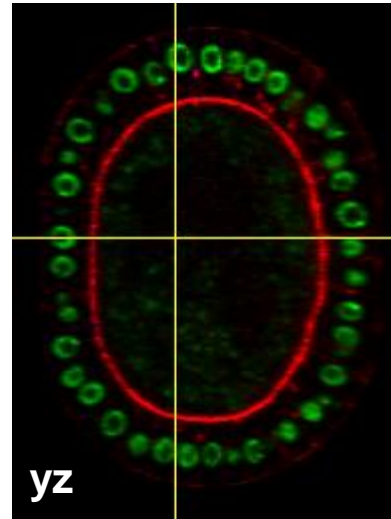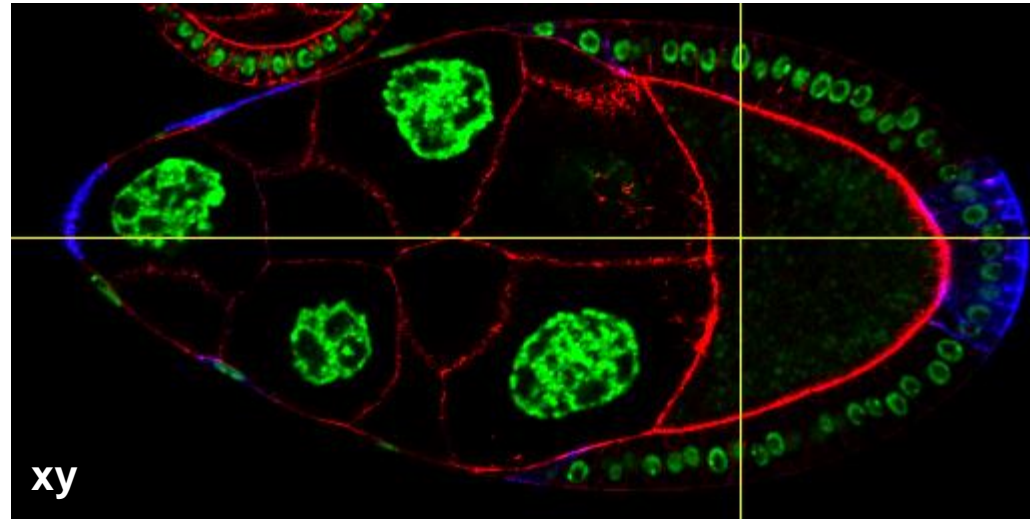


Orthogonal slice views

Example images

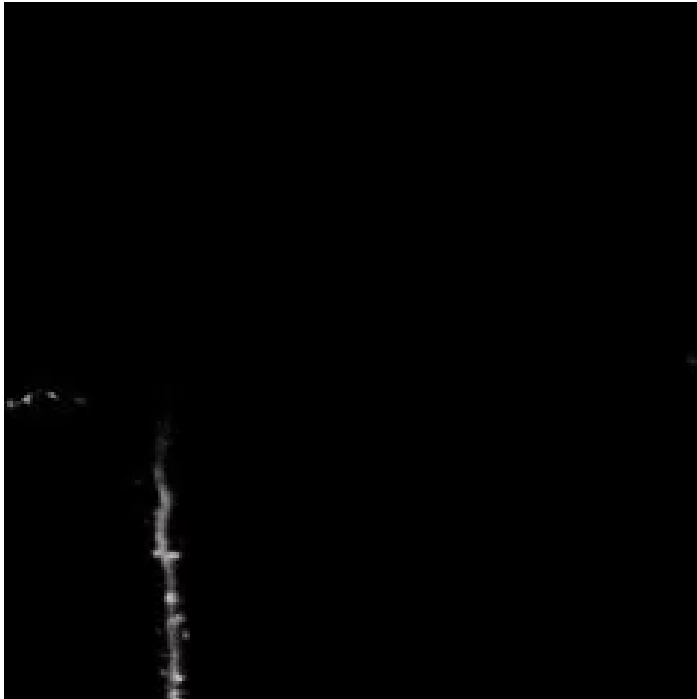

Imaris demo image

**Pyramidal cell z-stack
69 z-slices
15 frames-per-second playback**

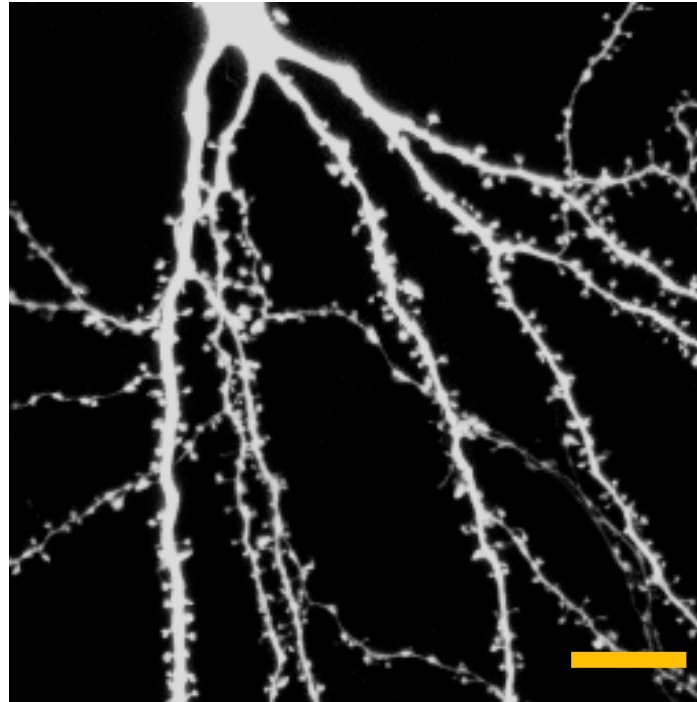


**Maximum intensity projection (Fiji)
Scalebar: 10 μ m**




3D reconstruction in Imaris


Example images

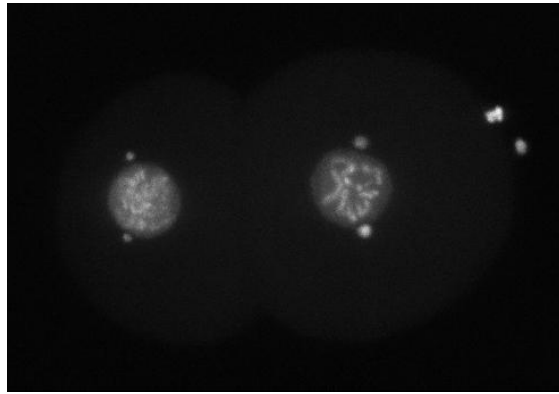


C. elegans embryo
eGFP-H2B
eGFP-gamma-tubulin
Max projection movie
One z-stack per minute acquisition
4 frames-per-second playback

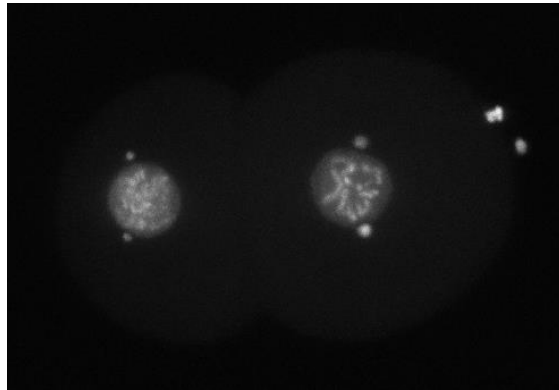


0 minutes

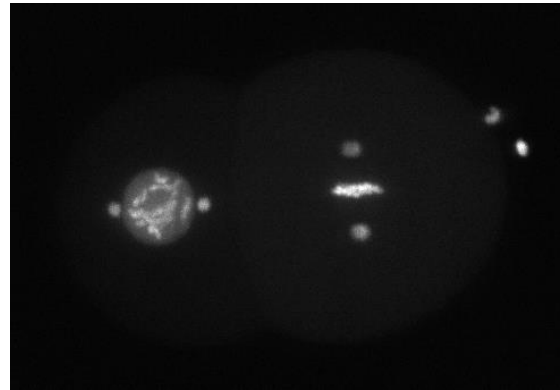

2 minutes

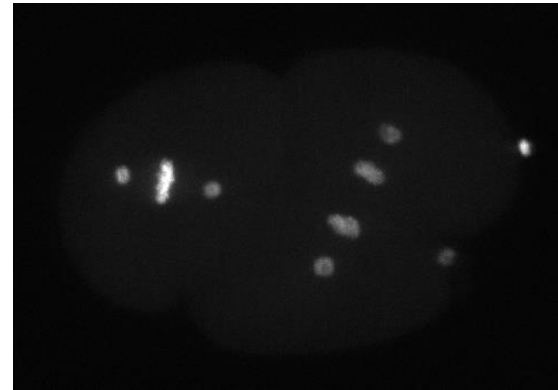


4 minutes

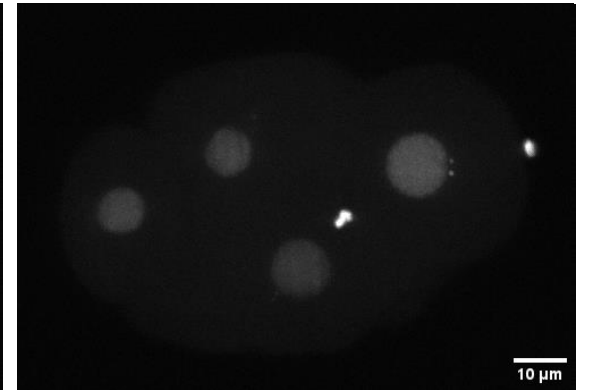


8 minutes

Example images

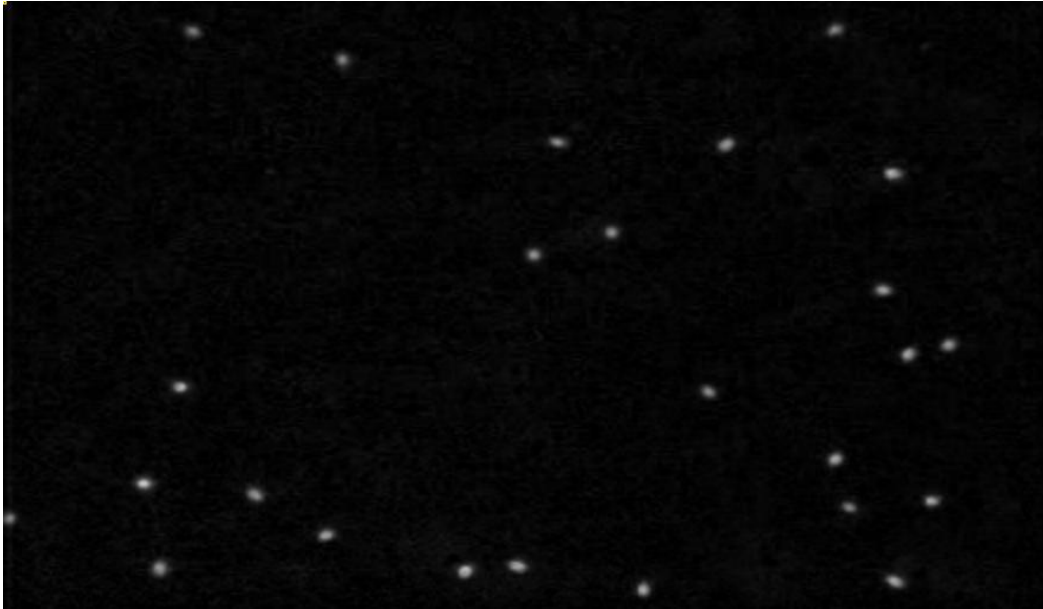


Imaris Demo Image – Swimming Algae
1 frame per second acquisition
12 frames per second playback

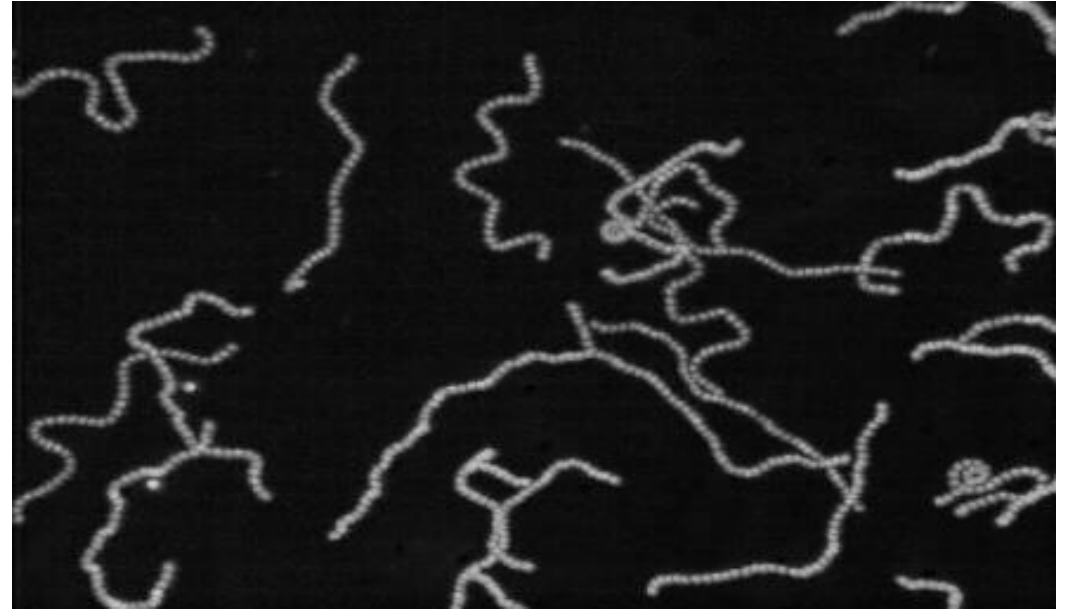


Temporal projection (Fiji)

Exercise 1 – raw vs composite RGB images

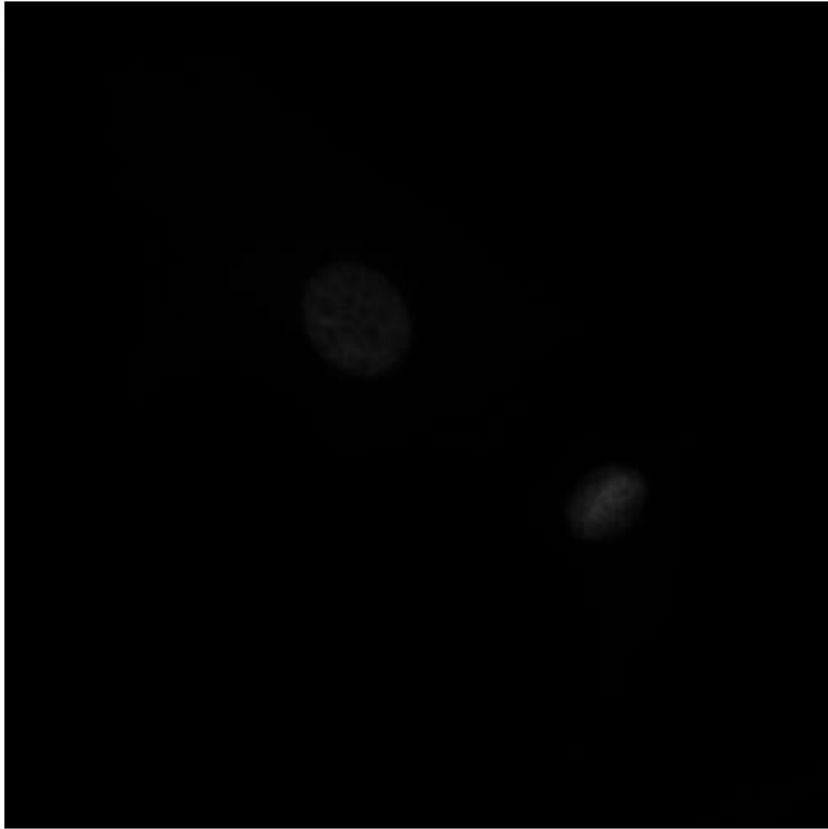

demo_DAPI_Phalloidin_Mitotracker_001.tif

Raw image, drag-and-drop in powerpoint.

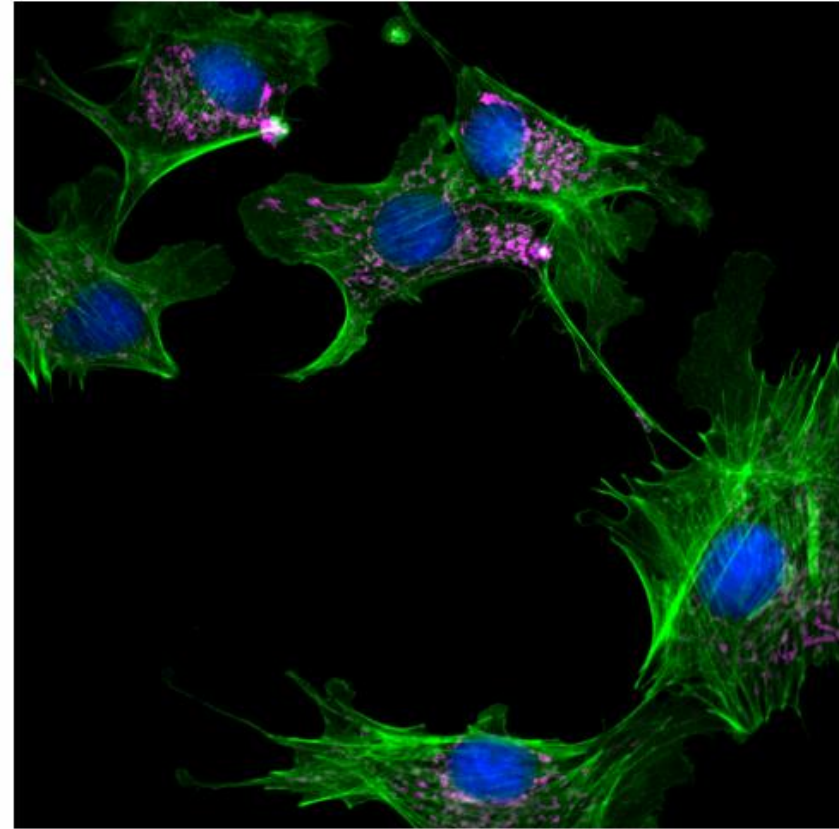

demo_DAPI_Phalloidin_Mitotracker_001.tif (RGB).tif


Composite RGB image generated with Exercise 1 macro,
drag-and-drop in powerpoint.

Day 2: Quantitative imaging analysis in Fiji


- Goal:
 - Learn to use common image processing tools in Fiji in order to extract quantitative information

Introduction to quantitative imaging analysis


- Why quantitative analysis?

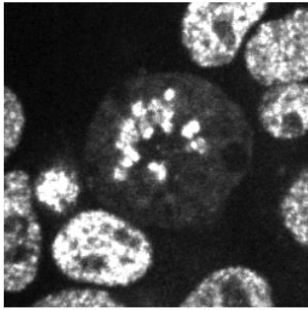

Tet1

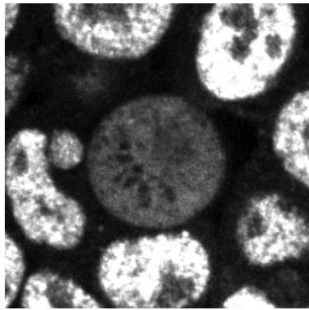


Tet1 $\Delta 1-131$

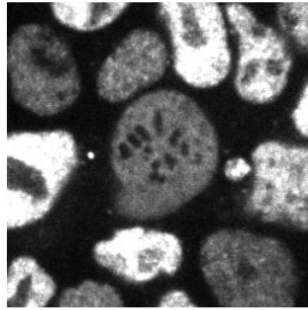


Tet2

Introduction to quantitative imaging analysis


- Why quantitative analysis?

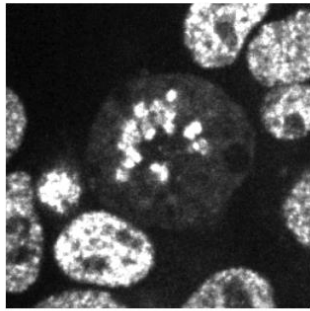


Tet1

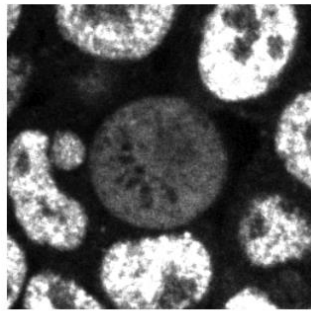

Tet1 $\Delta 1-131$

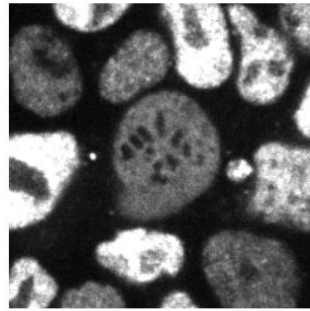


Tet2

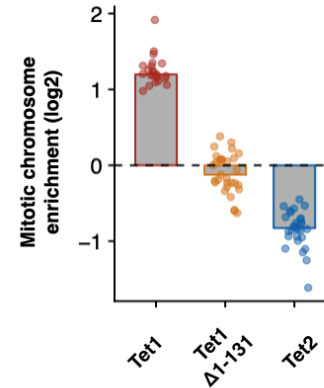


Introduction to quantitative imaging analysis


- Why quantitative analysis?

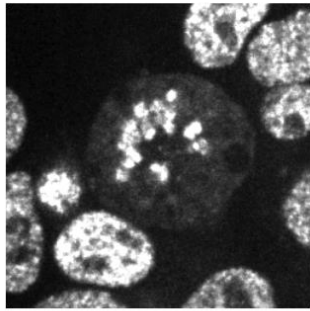


Tet1

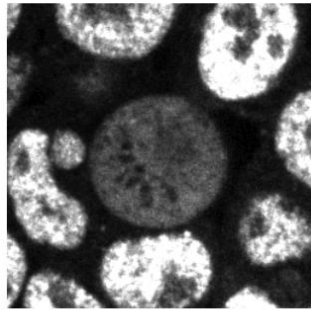

Tet1 Δ 1-131

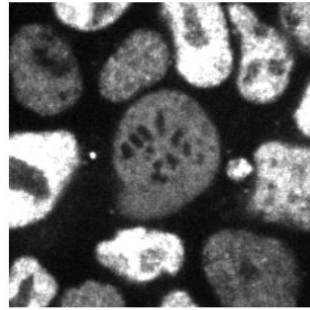

Tet2

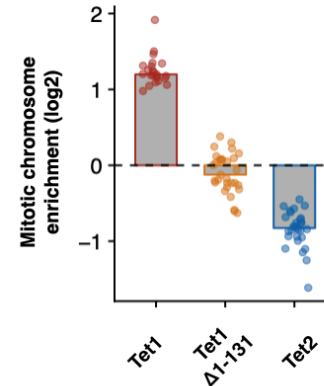

- Measured quantity (space, time, intensity etc)
- Increases robustness of observations
- Elucidate patterns that are not easily visible by eye




Better questions, better science!

Introduction to quantitative imaging analysis

Images can be quantitative:

- A digital image is simply a table of numbers
- Each number represents the intensity value of a pixel
- The intensity value is proportional to the amount of light collected at that location

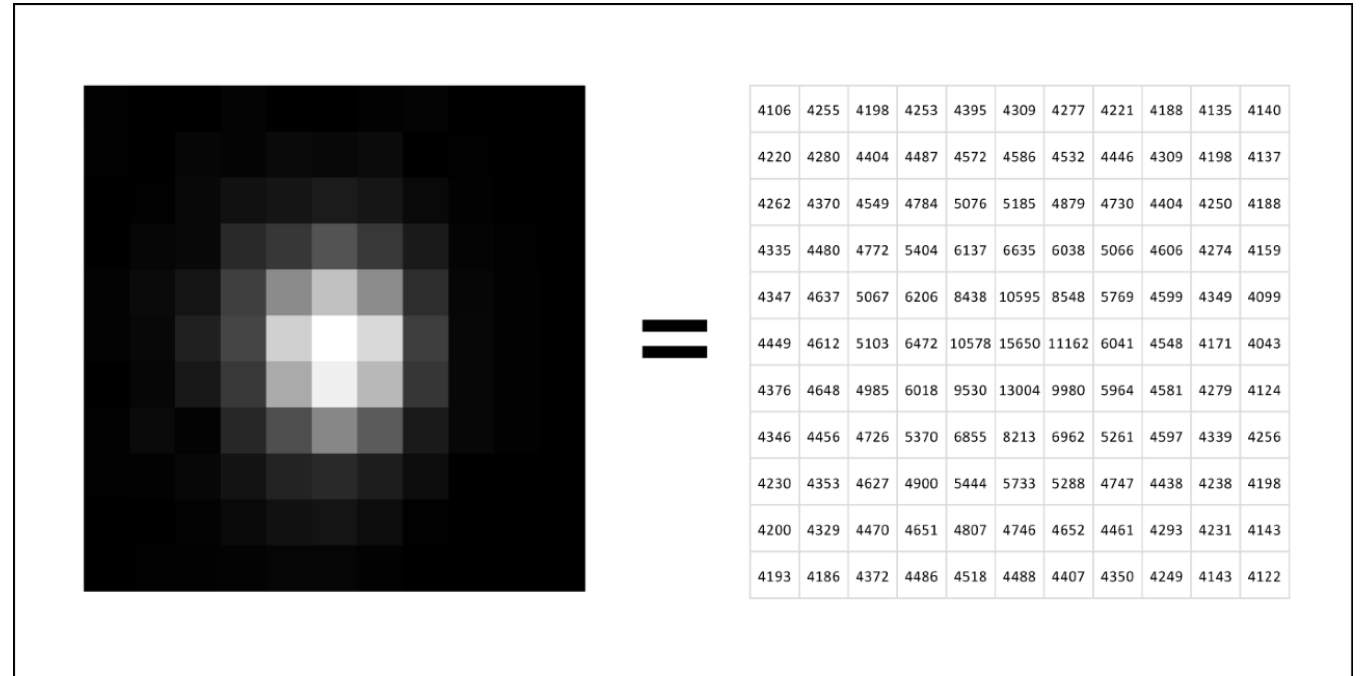


Value at each pixel = **Signal + Background + Noise**

Introduction to quantitative imaging analysis

- Goal:
 - Learn to use common image processing tools in Fiji in order to extract quantitative information

Introduction to quantitative imaging analysis


- Goal:
 - Learn to use common image processing tools in Fiji in order to extract quantitative information

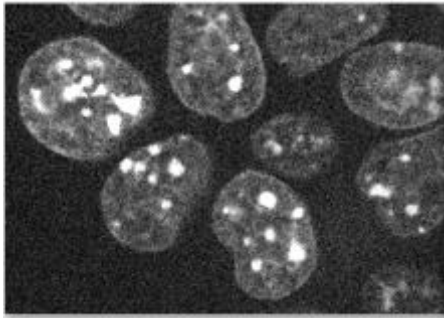


Open an image

Introduction to quantitative imaging analysis

- Goal:
 - Learn to use common image processing tools in Fiji in order to extract quantitative information

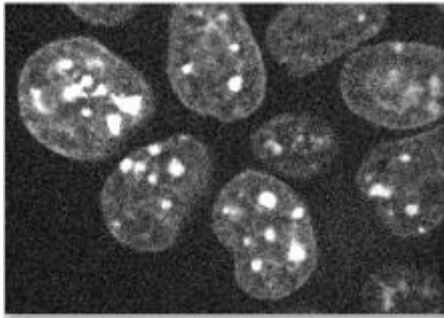


Open an image

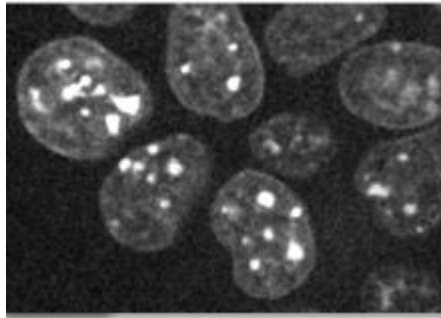


Processing such
as filtering

Introduction to quantitative imaging analysis

- Goal:
 - Learn to use common image processing tools in Fiji in order to extract quantitative information

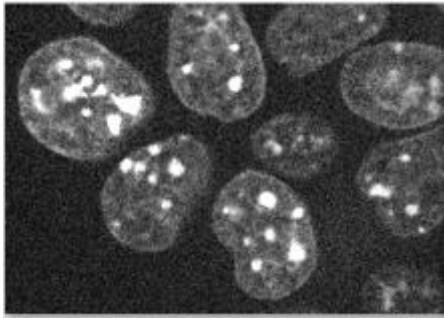

Open an image

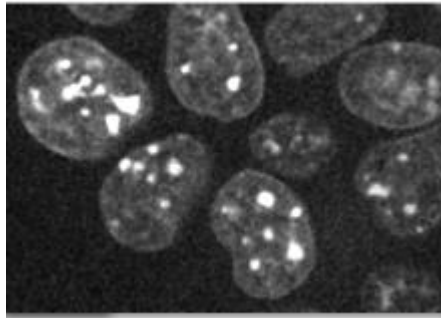


Processing such
as filtering




Segmentation or
object detection

Introduction to quantitative imaging analysis


- Goal:
 - Learn to use common image processing tools in Fiji in order to extract quantitative information

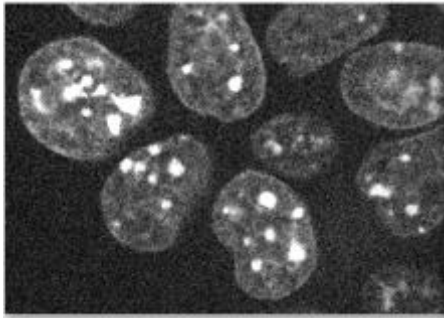


Open an image

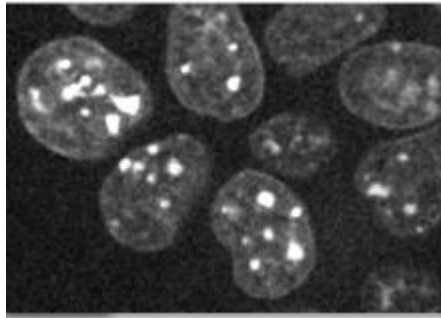

Processing such as filtering

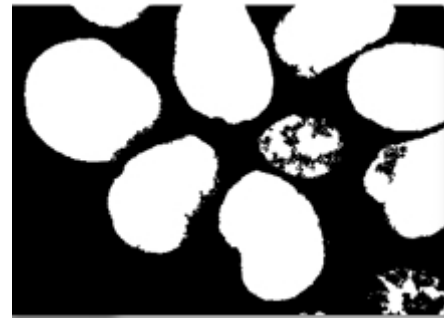

Segmentation or object detection

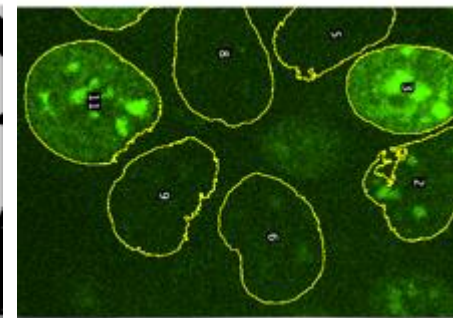

Data extraction from segmented objects

- Number of objects
- Shape / size
- Intensity within object...

Introduction to quantitative imaging analysis

- Goal:
 - Learn to use common image processing tools in **Fiji** in order to extract quantitative information




- Why Fiji?
 - Free open-source software
 - Huge development community
 - Libraries to open nearly any microscope image (!)
 - Automation and scripting with multiple programming languages (!)

Introduction to quantitative imaging analysis

- **Today:**

- Load images in Fiji, get familiarised with the program
- Processing using convolutional filters
- Segmentation with thresholding
- Data extraction
 - ...plotting?

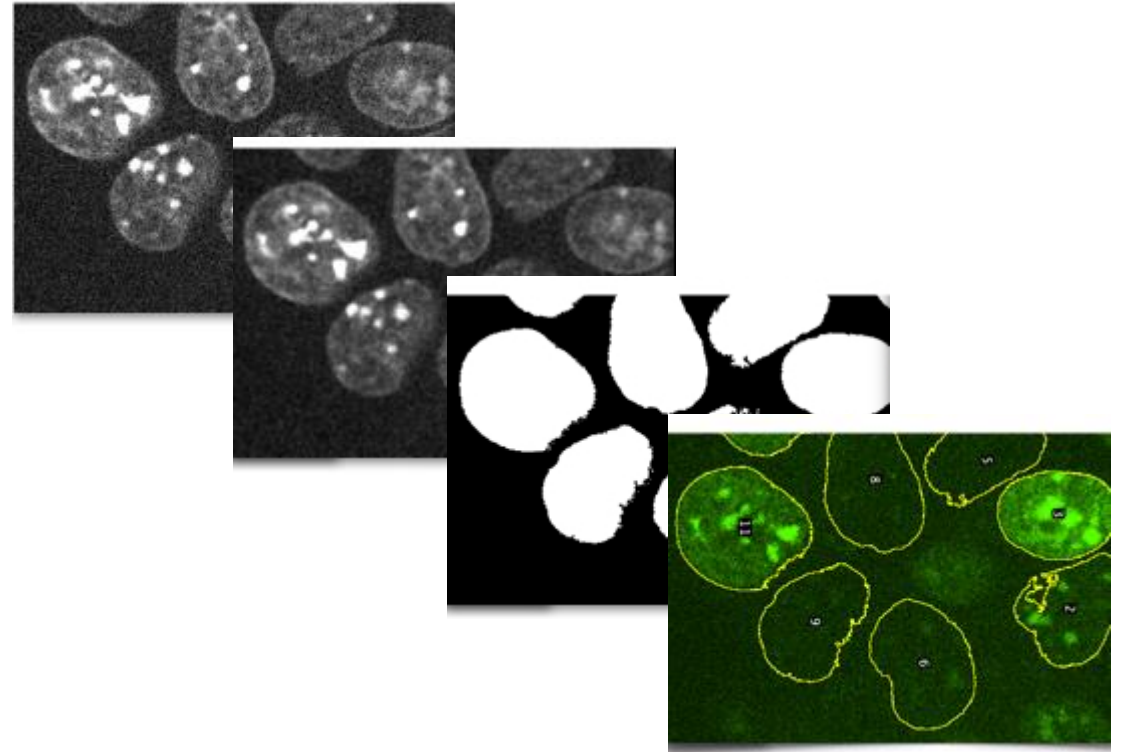

for loop

```
for (i = 0; i < numROI; i++) {  
    roiManager("Select", i);  
    Stack.setChannel(2);  
    run("Measure");  
}
```

for loop

Start value

End value

Increment

```
for (i = 0; i < numROI; i++) {  
  
    roiManager("Select", i);  
    Stack.setChannel(2);  
    run("Measure");  
  
}
```

for loop

Start value

End value

Increment

first, “build a sequence of integers going from the Start Value, to the End value, with an increment of 1”

```
for (i = 0; i < numROI; i++) {  
  
    roiManager("Select", i);  
    Stack.setChannel(2);  
    run("Measure");  
  
}
```

for loop

Start value

End value

Increment

```
for (i = 0; i < numROI; i++) {  
  
    roiManager("Select", i);  
    Stack.setChannel(2);  
    run("Measure");  
  
}
```

first, “build a sequence of integers going from the Start Value, to the End value, with an increment of 1”

---in the case of numROI being 5:

0, 1, 2, 3, 4

for loop

Start value

End value

Increment

```
for (i = 0; i < numROI; i++) {  
  
    roiManager("Select", i);  
    Stack.setChannel(2);  
    run("Measure");  
  
}
```

first, “build a sequence of integers going from the Start Value, to the End value, with an increment of 1”

---in the case of numROI being 5:

0, 1, 2, 3, 4,

“Run the { content of the loop } replacing i with the first integer of the sequence.

Then run the { content of the loop } replacing i with the second integer in the sequence... “ and so on...

for loop

Start value

End value

Increment

```
for (i = 0; i < numROI; i++) {
```

```
    roiManager("Select", i);
```

```
    Stack.setChannel(2);
```

```
    run("Measure");
```

```
}
```

```
i = 0  
roiManager("Select", i);  
Stack.setChannel(2);  
run("Measure");
```

```
i = 1  
roiManager("Select", i);  
Stack.setChannel(2);  
run("Measure");
```


```
i = 2  
roiManager("Select", i);  
Stack.setChannel(2);  
run("Measure");
```

```
i = 3  
roiManager("Select", i);  
Stack.setChannel(2);  
run("Measure");
```



```
...
```


Deep Learning based segmentation

- Sometimes objects are too close for adequate thresholding

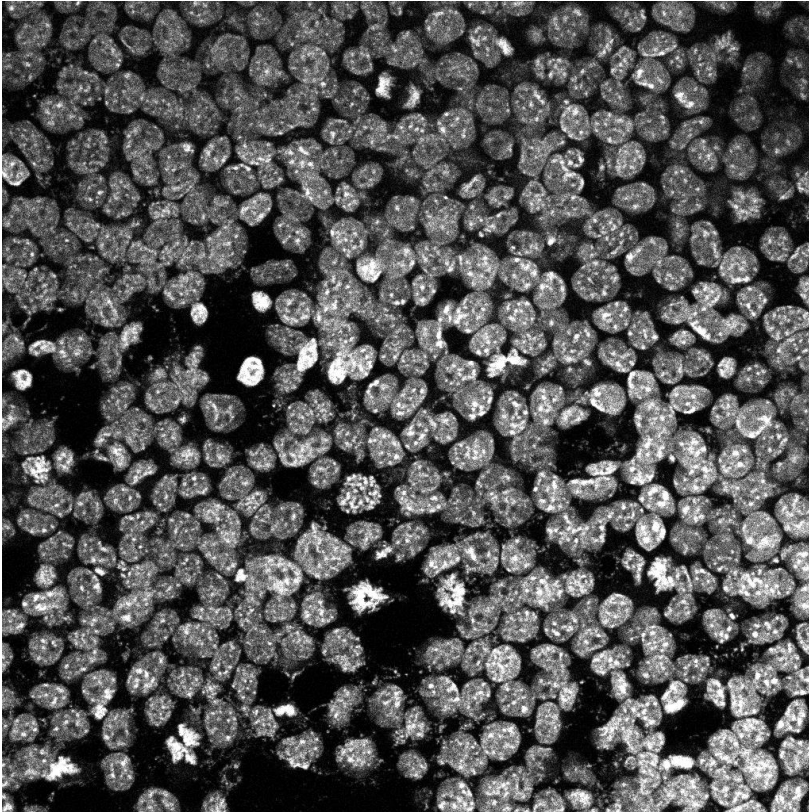


DAPI image

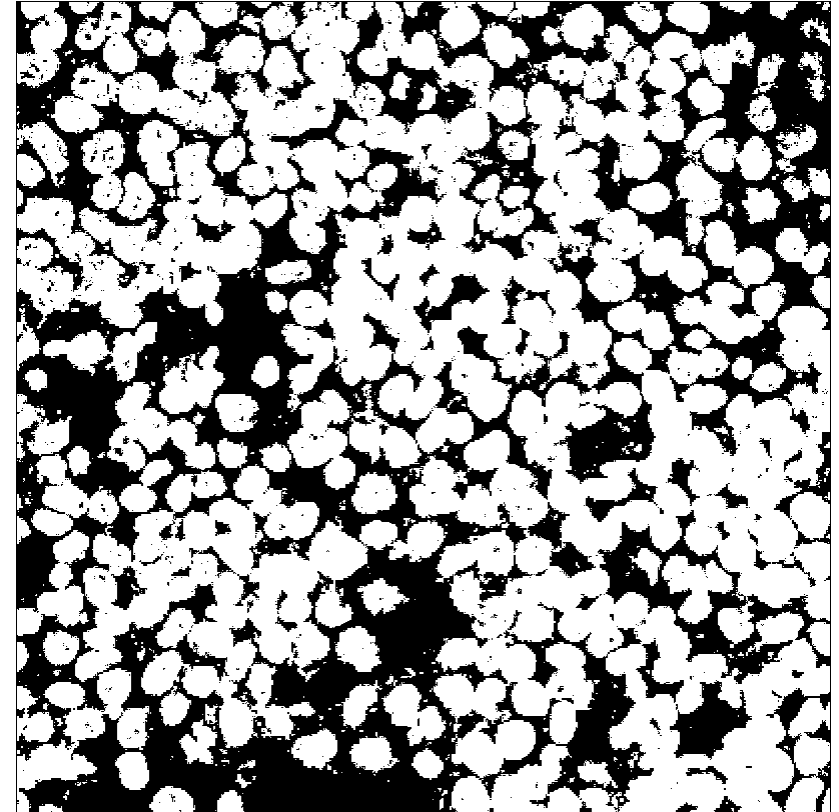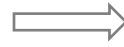

Otsu threshold image

Deep Learning based segmentation


- Machine learning tools trained on cell images can help distinguish close objects

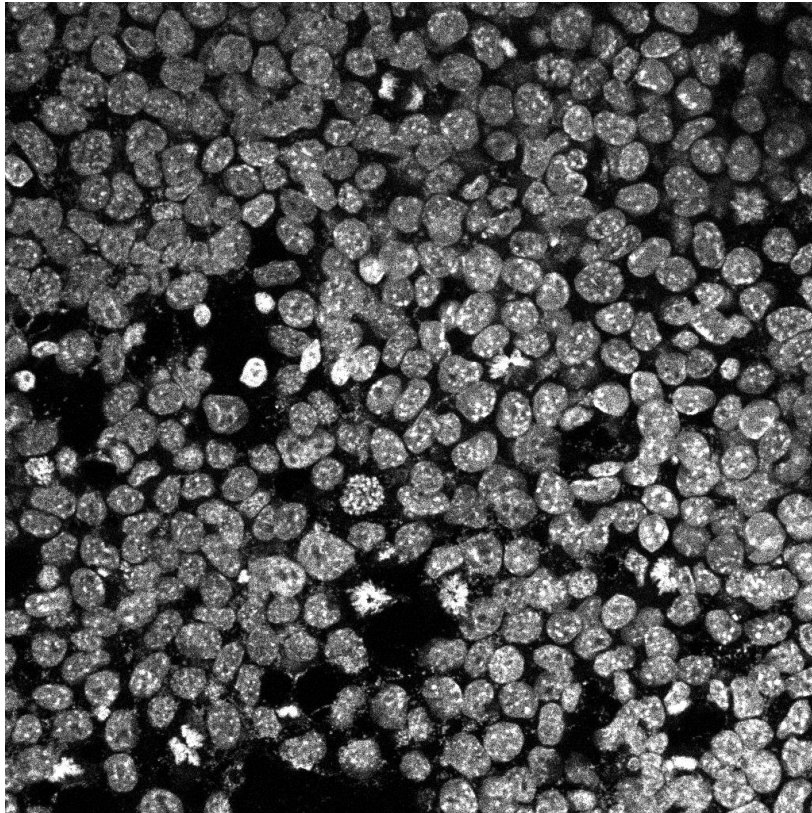

DAPI image




Nuclei masks identified with
Cellpose

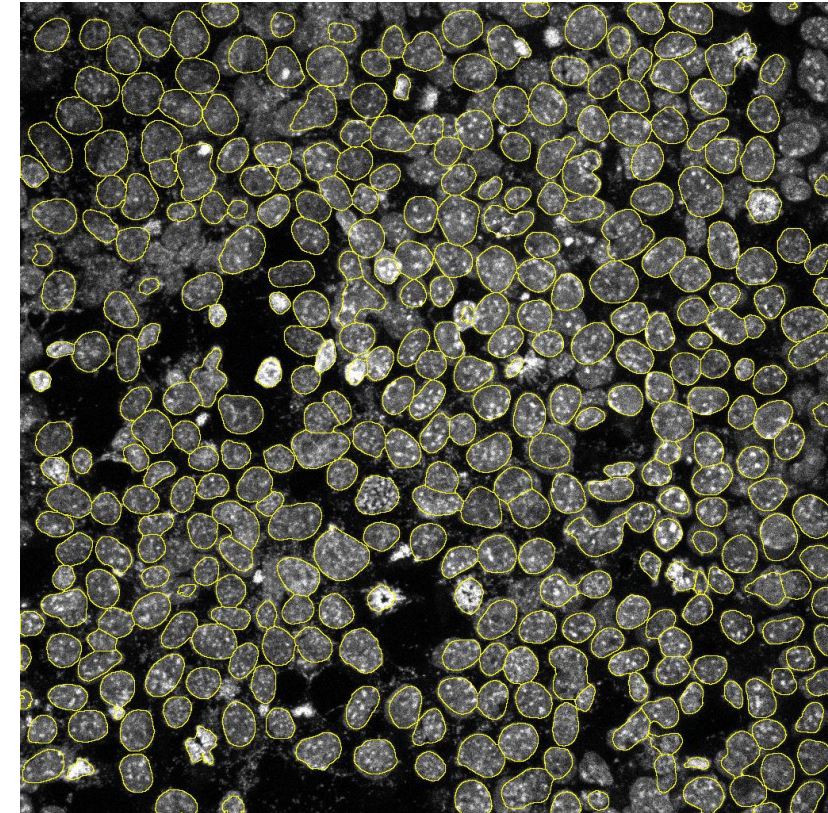


DAPI image + mask outlines

Deep Learning based segmentation


- Using the DAPI masks, we can then quantify intensities of other stainings in the nucleus

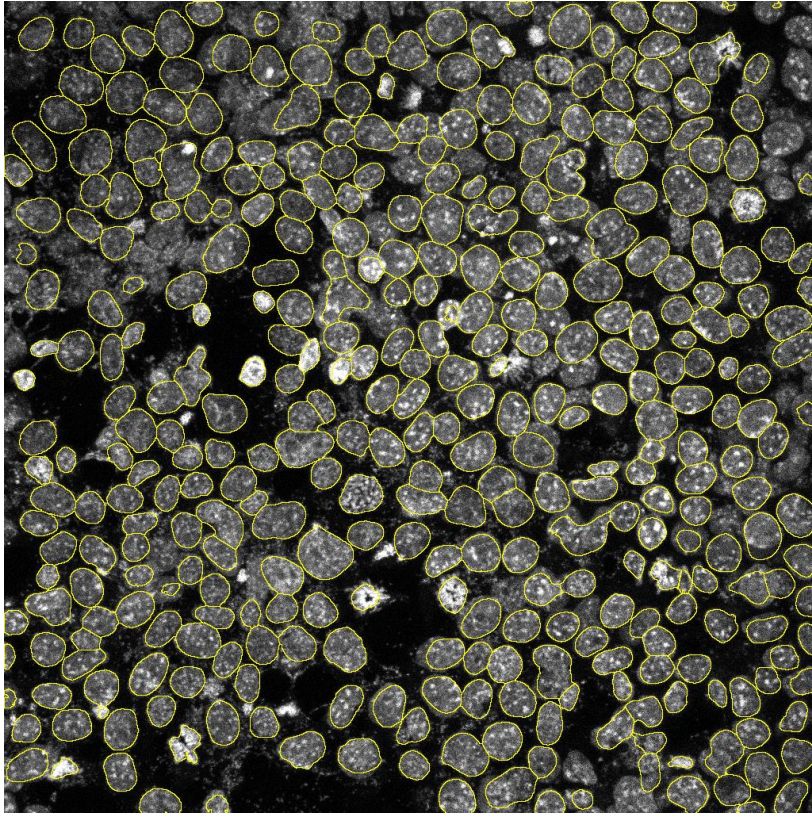

DAPI image + nucleus mask outlines

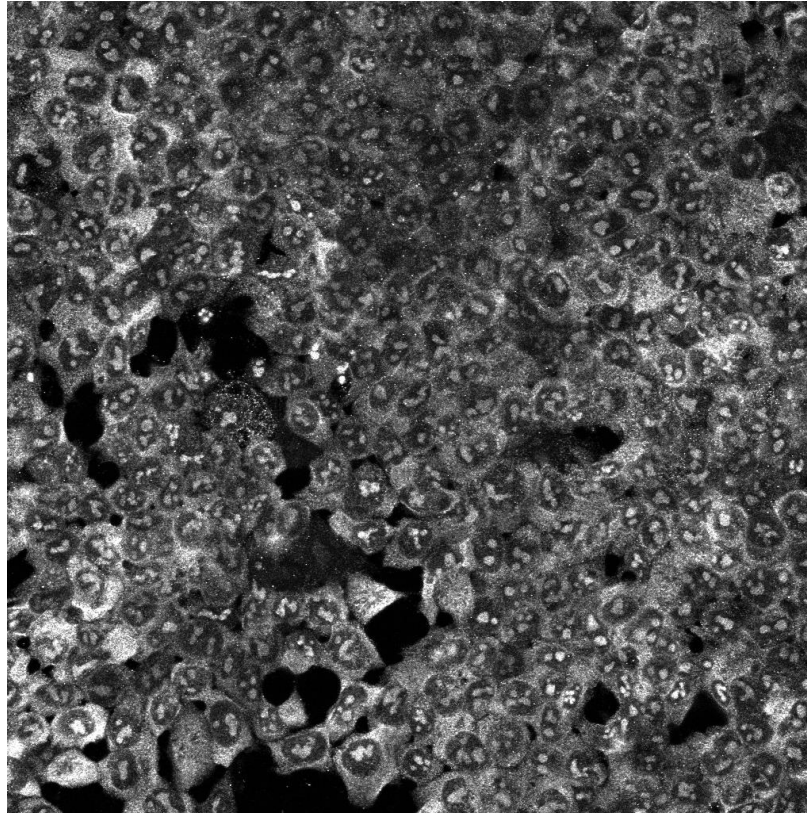


LIN28-GFP image

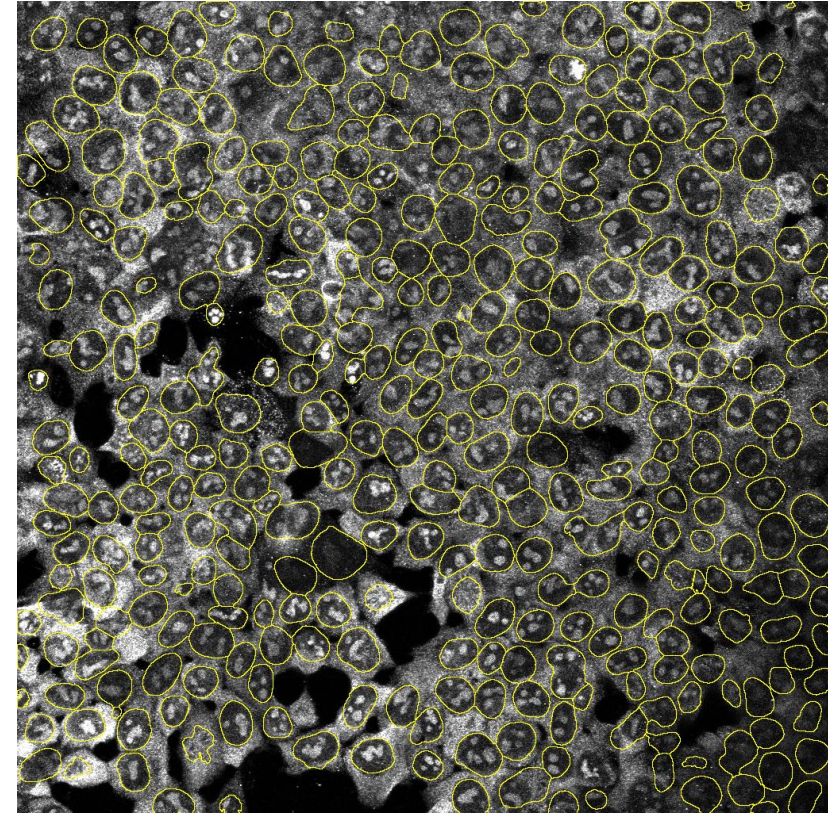

LIN28-GFP image + nucleus mask outlines


Deep Learning based segmentation

- Using the DAPI masks, we can then quantify intensities of other stainings in the nucleus

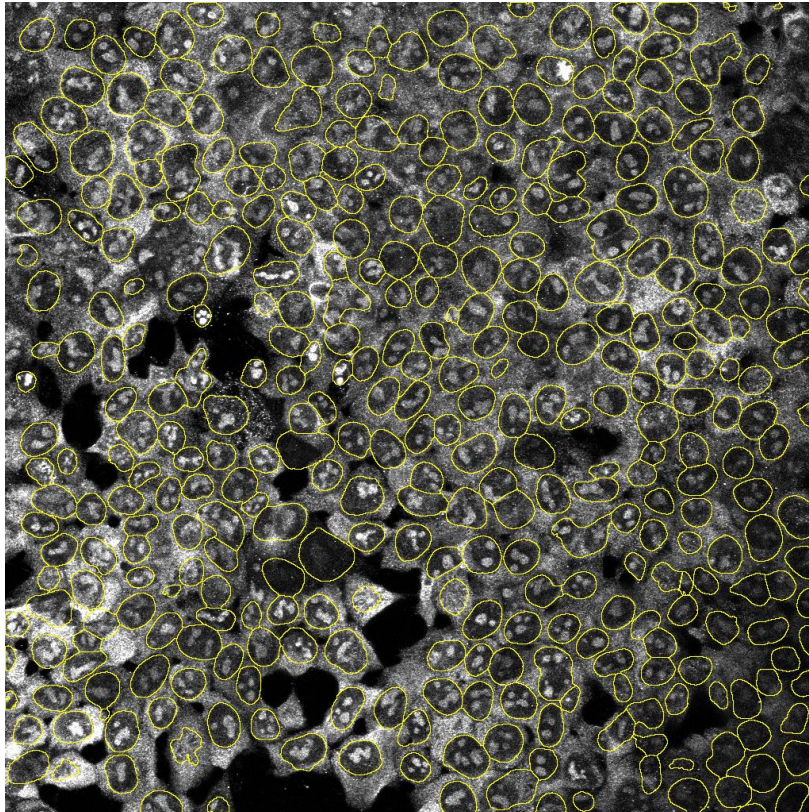

LIN28-GFP image + nucleus mask outlines

Results									
File Edit Font Results									
Label	Area	Mean	StdDev	Mode					
342 MAX_20240129_CGR8_iBidi_Timecourse_OTX2_SOX2_ESRRB_LIN28A_CGR8_OTX2_T0_1-2.tif 0949-0908	147.453	12.157	5.874	10					
343 MAX_20240129_CGR8_iBidi_Timecourse_OTX2_SOX2_ESRRB_LIN28A_CGR8_OTX2_T0_1-2.tif 0944-0376	155.523	32.527	17.836	20					
344 MAX_20240129_CGR8_iBidi_Timecourse_OTX2_SOX2_ESRRB_LIN28A_CGR8_OTX2_T0_1-2.tif 0935-0813	66.584	12.790	7.279	10					
345 MAX_20240129_CGR8_iBidi_Timecourse_OTX2_SOX2_ESRRB_LIN28A_CGR8_OTX2_T0_1-2.tif 0950-0581	187.322	26.753	14.969	20					
346 MAX_20240129_CGR8_iBidi_Timecourse_OTX2_SOX2_ESRRB_LIN28A_CGR8_OTX2_T0_1-2.tif 0951-0448	148.421	36.044	22.070	20					
347 MAX_20240129_CGR8_iBidi_Timecourse_OTX2_SOX2_ESRRB_LIN28A_CGR8_OTX2_T0_1-2.tif 0952-1002	116.864	8.281	4.259	5					
348 MAX_20240129_CGR8_iBidi_Timecourse_OTX2_SOX2_ESRRB_LIN28A_CGR8_OTX2_T0_1-2.tif 0960-0678	87.810	22.566	14.736	11					
349 MAX_20240129_CGR8_iBidi_Timecourse_OTX2_SOX2_ESRRB_LIN28A_CGR8_OTX2_T0_1-2.tif 0963-0029	133.087	27.723	14.907	15					
350 MAX_20240129_CGR8_iBidi_Timecourse_OTX2_SOX2_ESRRB_LIN28A_CGR8_OTX2_T0_1-2.tif 0963-0517	83.694	29.380	21.045	12					
351 MAX_20240129_CGR8_iBidi_Timecourse_OTX2_SOX2_ESRRB_LIN28A_CGR8_OTX2_T0_1-2.tif 0971-0955	132.118	7.816	3.847	5					
352 MAX_20240129_CGR8_iBidi_Timecourse_OTX2_SOX2_ESRRB_LIN28A_CGR8_OTX2_T0_1-2.tif 0972-0146	127.518	30.303	17.061	20					
353 MAX_20240129_CGR8_iBidi_Timecourse_OTX2_SOX2_ESRRB_LIN28A_CGR8_OTX2_T0_1-2.tif 0972-0813	83.774	16.490	11.709	7					
354 MAX_20240129_CGR8_iBidi_Timecourse_OTX2_SOX2_ESRRB_LIN28A_CGR8_OTX2_T0_1-2.tif 0959-0852	25.905	13.484	0.645	8					
355 MAX_20240129_CGR8_iBidi_Timecourse_OTX2_SOX2_ESRRB_LIN28A_CGR8_OTX2_T0_1-2.tif 0967-0076	44.793	60.728	21.658	66					
356 MAX_20240129_CGR8_iBidi_Timecourse_OTX2_SOX2_ESRRB_LIN28A_CGR8_OTX2_T0_1-2.tif 0991-0640	192.165	28.741	13.635	19					
357 MAX_20240129_CGR8_iBidi_Timecourse_OTX2_SOX2_ESRRB_LIN28A_CGR8_OTX2_T0_1-2.tif 0994-0918	115.089	12.575	5.527	8					

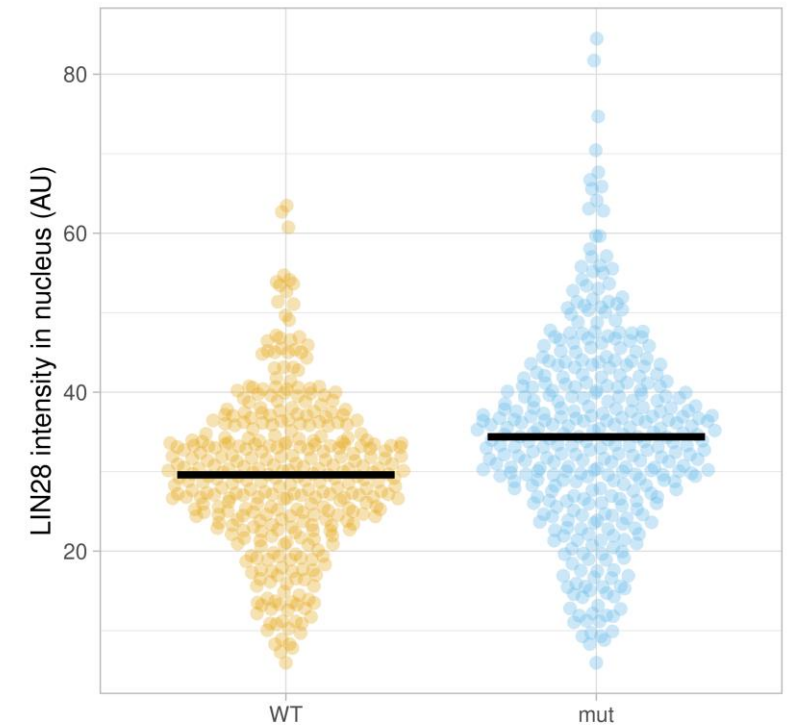


Plotting resources

PlotsOfData - Plots all Of the Data

Data upload

- ☒ Example 1 (wide format)
- ☐ Example 2 (tidy format)
- ☐ Upload file
- ☐ Paste data
- ☐ URL (csv files only)
- ☐ These data are Tidy

Select and order:

 Download in tidy format (csv)

☐ Show information on data formats

Data offset

- ☒ Quasirandom
- ☐ Random
- ☐ None; stripes
- ☐ None (for small n)

Visibility of the data


0

0.3

1

Statistics

- ☒ Median
- ☐ Mean
- ☐ Boxplot (minimum n=10)
- ☐ Violin Plot (minimum n=10)
- ☐ Add 95% CI (minimum n=10)

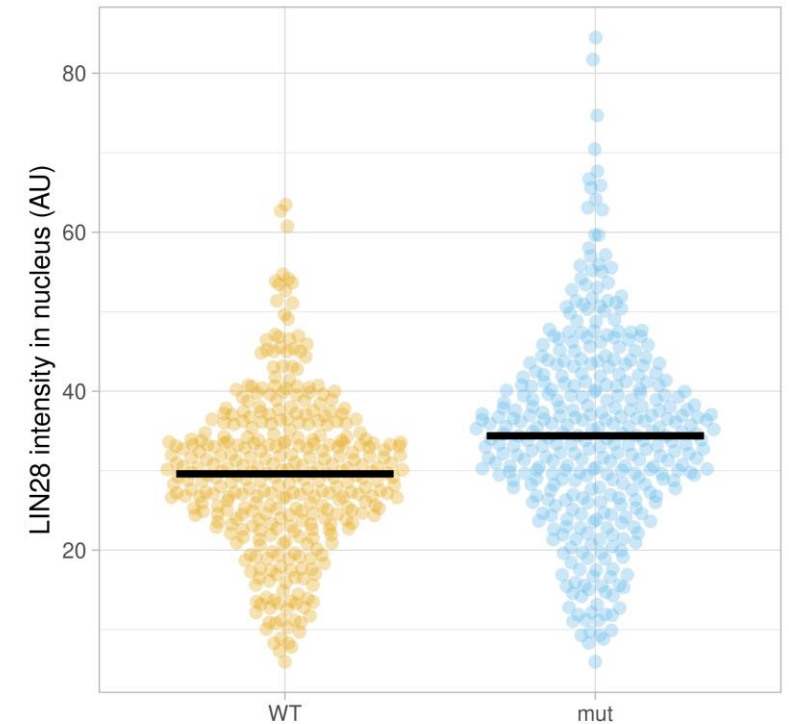

Generate plots with data points on PlotsOfData webapp

<https://huygens.science.uva.nl/PlotsOfData/>